

easydrip™ Plus

Pen Needles

Innovation for easy & conformable injections

New generation pen needles provide with more comfort against the skin during the injection process.



Reliable Injection Depth



High Flow Rate




Silicone Coating



New Generation Design

More comfortable against the skin during the injection process

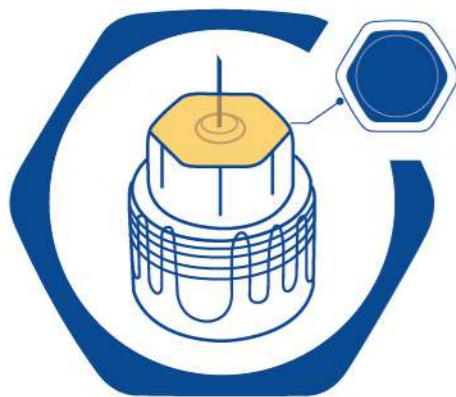


*Pen injectors compatible with the pen needle: NovoPen®5, Autopen®24, Gansulin Pen, Ypsopen®, UNIPEN®, XiuLin Pen, HumaPen®ErgoII, Byetta, Victoza Pen
* All the brands are trademarks of their respective owners. Please read insulin pen instructions for use.

Contoured hub for Reliable Depth Control

The contoured hub distributes insertion forces across a wider skin surface with the intent to minimize the clinical impact (IM injection) of variable patient injection forces on injection depot depth.


Intramuscular (IM) injection can increase insulin absorption, causing hypoglycemia.

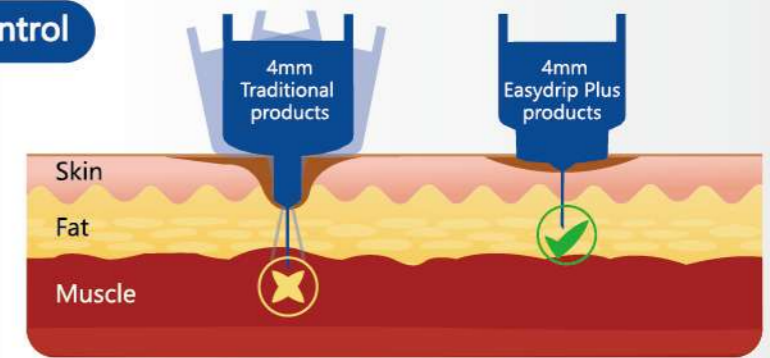


Ultra-strong Design to resist bending or breaking



An improved ability to hold the pen securely against the skin without wobbling during an injection.

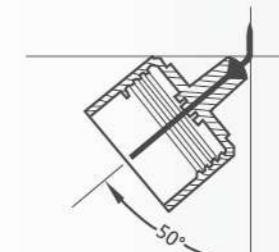

Traditional products

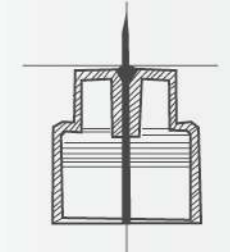

Easydrip Plus products

Ergonomic design helping diabetes use PN's more easily

Wider outer cover with ergonomically designed grooves and a colored, larger, more prominent inner shield makes it easier to remove the inner shield.



Wider inner cap



Wider outer cap



Ordering information

Reorder #	Color	Needle gauge	Needle length	Standard packing
02-3204		32G needle	4 mm	100 pcs/box 2,000 pcs/case
02-3105		31G needle	5 mm	
02-3106		31G needle	6 mm	

1. Needle gauge refers to the size of the needle. The higher the gauge, the thinner the needle;
2. Please contact our sales person if you have special requests for an unique size.